

# MAGLUMI LH (CLIA)



130202002M



100



**Shenzhen New Industries  
Biomedical Engineering Co., Ltd**  
4/F, Wearnes Tech Bldg,  
Science & Industry Park,  
Nanshan, Shenzhen, 518057 CHINA  
Tel. + 86-755-86028224  
Fax. + 86-755-26654850



**Lotus Global Co., Ltd**  
15 Alexandra Road  
London  
UK  
NW8 0DP  
Tel. + 44-20-75868010  
Fax. + 44-20-79006187



## FOR PROFESSIONAL USE ONLY

Store at 2-8 °C



COMPLETELY READ THE INSTRUCTIONS BEFORE  
PROCEEDING



## SYMBOLS EXPLANATIONS



Authorized Representative in the  
European community



Manufacturer



Consult instructions for use



Contents of kit



In vitro diagnostic medical device




Batch code



Catalogue number



Use by

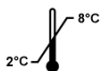

Temperature limitation  
(store at 2-8 °C)



Sufficient for



Keep away from sunlight



Keep upright for storage

## INTENDED USE

The kit has been designed for the quantitative determination of luteal hormone (LH) in human serum.

The method can be used for samples over the range of 0.5-250 mIU/ml.

The test has to be performed on the Fully-auto chemiluminescence immunoassay (CLIA) analyzer MAGLUMI (Including Maglumi 600, Maglumi 1000, Maglumi 1000 Plus, Maglumi 2000, Maglumi 2000 Plus, Maglumi 3000 and Maglumi 4000).

## SUMMARY AND EXPLANATION OF THE TEST

The gonadotrophins LH and FSH glycoproteins with a molecular weight of approx. 30,000 Dalton are secreted by the basophilic cells of the pituitary gland. Their pulsatile secretion is regulated by the hypothalamic gonadotrophin releasing hormone (GnRH, LHRH).

In the female, the gonadotrophins stimulate the growth of ovarian follicles during the follicular phase. During the luteal phase, their secretion is inhibited by the influence of progesterone and estradiol. In menopausal women, there is a striking rise in FSH serum concentrations, while LH does not exceed the usual level of the mid-cycle peak. High postmenopausal LH and FSH levels are due to the absence of the progesterone/ estradiol response.

In the male, LH stimulates the testosterone production of the Leydig cells. In combination with FSH and testosterone, it also stimulates spermatogenesis. The secretion of LH is regulated by a negative feedback mechanism of testosterone; the regulation of FSH is subject to the influence of inhibin.

During childhood, LH and FSH levels are normally too low to be detectable. With the beginning of puberty, FSH is the first gonadotrophin to reach detectable values.

The determination of LH and FSH plays an important role in the detection of dysfunctions of the pituitary-ovarian axis, clinically manifested by amenorrhoea, oligomenorrhoea, anovulatory cycles and menorrhagia. For differential diagnosis between hypothalamic and pituitary disturbances, the GnRH and the clomiphene tests are widely used. Determination of LH serum concentrations is also used as an indicator of ovulation in in-vitro fertilisation.

In men, determination of gonadotrophins is mainly used for the differentiation between primary and secondary hypogonadism.

In children and juveniles, LH and FSH determinations are indicated in case of either precocious or delayed puberty.

## PRINCIPLE OF THE TEST

Sandwich immunoluminometric assay:

Use an anti-LH monoclonal antibody to label ABEI, and use another monoclonal antibody to coat magnetic microbeads. Sample, Calibrator or Control, with ABEI Label and magnetic microbeads are mixed thoroughly and incubated at 37°C, forming a sandwich; After sediment in a magnetic field, suck the supernatant and then cycle wash it for 1 time. Subsequently, the starter reagents are added and a flash chemiluminescent reaction is initiated. The light signal is measured by a photomultiplier as RLU within 3 seconds and is proportional to the concentration of LH present in samples.



## KIT COMPONENTS

### Material Supplies

Reagent Integral for 100 determinations	
Nano magnetic microbeads: TRIS buffer, 1.2%(W/V), 0.2%NaN <sub>3</sub> , coated with anti-LH monoclonal antibody.	2.5ml

<b>Calibrator Low:</b> bovine serum, 0.2%NaN <sub>3</sub>	2.5ml
<b>Calibrator High:</b> bovine serum, 0.2%NaN <sub>3</sub>	2.5ml
<b>ABEI Label:</b> anti-LH monoclonal antibody labeled ABEI, containing BSA, 0.2%NaN <sub>3</sub> .	10.5ml
All reagents are provided ready-to-use.	

Reagent Vials in kit box	
<b>Internal Quality Control:</b> containing BSA, 0.2%NaN <sub>3</sub> . (target value refer to Quality Control Information date sheet)	2.0ml

Internal quality control is only applicable with MAGLUMI system. Instructions for use and target value refer to Quality Control Information date sheet. User needs to judge results with their own standards and knowledge.

#### Accessories Required But Not Provided

MAGLUMI Reaction Module	REF: 630003
MAGLUMI Starter 1+2	REF: 130299004M
MAGLUMI Wash Concentrate	REF: 130299005M
MAGLUMI Light Check	REF: 130299006M

Please order accessories from SNIBE or our representative.



#### Preparation of the Reagent Integral

Before the sealing is removed, gentle and careful horizontal shaking of the Reagent Integral is essential (avoid foam formation!) Remove the sealing and turn the small wheel of the magnetic microbeads compartment to and fro, until the colour of the suspension has changed into brown. Place the Integral into the reagent area and let it stand there for 30 min. During this time, the magnetic microbeads are automatically agitated and completely resuspended.

**Do not interchange integral component from different reagents or lots!**

#### Storage and Stability

- Sealed: Stored at 2-8°C until the expiry date.
- Opened: Stable for 4 weeks. To ensure the best kit performance, it is recommended to place opened kits in the refrigerator if it's not going to be used on board during the next 12 hours.



- Keep upright for storage.



- Keep away from sunlight.

## CALIBRATION AND TRACEABILITY

### 1) Traceability

To perform an accurate calibration, we have provided the test calibrators standardized against the W.H.O 2nd International Standard 80/552

### 2) 2-Point Recalibration

Via the measurement of calibrators, the predefined master curve is adjusted (recalibrated) to a new, instrument-specific measurement level with each calibration.

### 3) Frequency of Recalibration

- After each exchange of lot (Reagent Integral or Starter Reagents).
  - Every 4 weeks and/or each time a new Integral is used (recommendation).
  - After each servicing of the Fully-auto chemiluminescence immunoassay (CLIA) analyzer MAGLUMI.
  - If controls are beyond the expected range.
  - The room temperature has changed more than 5 °C
- 012130909-V2.3-EN

(recommendation).

## SPECIMEN COLLECTION AND PREPARATION

Sample material: serum

Collect 5.0ml venous blood into Blood Collection Tube. Standing at room temperature, centrifuging, separating serum part.

The serum sample is stable for up to 12 hours at 2-8°C. More than 12 hours, please packed, -20 °C can be stored for 30 days.

Avoid repeated freezing and thawing, the serum sample can be only frozen and thawed two times. Stored samples should be thoroughly mixed prior to use (Vortex mixer).

Please ask local representative of SNIBE for more details if you have any doubt.

#### Vacuum Tubes

- Blank tubes are recommended type for collecting samples.
- Please ask SNIBE for advice if special additive must be used in sample collecting.

#### Specimen Conditions

- Do not use specimens with the following conditions:
  - heat-inactivated specimens;
  - Cadaver specimens or body fluids other than human serum;
  - Obvious microbial contamination.
- Use caution when handling patient specimens to prevent cross contamination. Use of disposable pipettes or pipette tips is recommended.
- Inspect all samples for bubbles. Remove bubbles with an applicator stick prior to analysis. Use a new applicator stick for each sample to prevent cross contamination.
- Serum specimens should be free of fibrin, red blood cells or other particulate matter.
- Ensure that complete clot formation in serum specimens has taken place prior to centrifugation. Some specimens, especially those from patients receiving anticoagulant or thrombolytic therapy, may exhibit increased clotting time. If the specimen is centrifuged before a complete clot forms, the presence of fibrin may cause erroneous results.

#### Preparation for Analysis

- Patient specimens with a cloudy or turbid appearance must be centrifuged prior to testing. Following centrifugation, avoid the lipid layer (if present) when pipetting the specimen into a sample cup or secondary tube.
- Specimens must be mixed **thoroughly** after thawing by **low** speed vortexing or by gently inverting, and centrifuged prior to use to remove red blood cells or particulate matter to ensure consistency in the results. Multiple freeze-thaw cycles of specimens should be avoided.
- All samples (patient specimens or controls) should be tested within 3 hours of being placed on board the MAGLUMI System. Refer to the SNIBE service for a more detailed discussion of onboard sample storage constraints.

#### Storage


- If testing will be delayed for more than 8 hours, remove serum or plasma from the serum separator, red blood cells or clot. Specimens removed from the separator gel, cells or clot may be stored up to 12 hours at 2-8°C.
- Specimens can be stored up to 30 days frozen at -20°C or colder.

#### Shipping

- Before shipping specimens, it is recommended that specimens be removed from the serum separator, red blood cells or clot. When shipped, specimens must be packaged and labeled in compliance with applicable state, federal and international regulations covering the transport of clinical specimens and infectious substances. Specimens must be shipped frozen

(dry ice). Do not exceed the storage time limitations identified in this section of the package insert.

## WARNING AND PRECAUTIONS FOR USERS

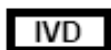

- For use in *IN-VITRO* diagnostic procedures only.
- Package insert instructions must be carefully followed. Reliability of assay results cannot be guaranteed if there are any deviations from the instructions in this package insert.

### Safety Precautions

**CAUTION:** This product requires the handling of human specimens.

- The calibrators in this kit are prepared from bovine serum products. However, because no test method can offer complete assurance that HIV, Hepatitis B Virus or other infectious agents are absent; these reagents should be considered a potential biohazard and handled with the same precautions as applied to any serum or plasma specimen.
- All samples, biological reagents and materials used in the assay must be considered potentially able to transmit infectious agents. They should therefore be disposed of in accordance with the prevailing regulations and guidelines of the agencies holding jurisdiction over the laboratory, and the regulations of each country. Disposable materials must be incinerated; liquid waste must be decontaminated with sodium hypochlorite at a final concentration of 5% for at least half an hour. Any materials to be reused must be autoclaved using an overkill approach. A minimum of one hour at 121°C is usually considered adequate, though the users must check the effectiveness of their decontamination cycle by initially validating it and routinely using biological indicators.
- It is recommended that all human sourced materials be considered potentially infectious and handled in accordance with the OSHA Standard on Bloodborne Pathogens<sup>13</sup>. Biosafety Level 214 or other appropriate biosafety practices should be used for materials that contain or are suspected of containing infectious agents.
- This product contains Sodium Azide; this material and its container must be disposed of in a safe way.
- Safety data sheets are available on request.

### Handling Precautions

- Do not use reagent kits beyond the expiration date.
- Do not mix reagents from different reagent kits.
- Prior to loading the Reagent Kit on the system for the first time, the microbeads requires mixing to re-suspend microbeads that have settled during shipment.
- For microbeads mixing instructions, refer to the KIT COMPONENTS, Preparation of the Reagent Integral section of this package insert.
- To avoid contamination, wear clean gloves when operating with a reagent kit and sample.
- Over time, residual liquids may dry on the kit surface, please pay attention the silicon film still exists on the surface of the kit.
- For a detailed discussion of handling precautions during system operation, refer to the SNIBE service information.

## TEST PROCEDURE

To ensure proper test performance, strictly adhere to the operating instructions of the Fully-auto chemiluminescence immunoassay (CLIA) analyzer MAGLUMI. Each test parameter is identified via a RFID tag on the Reagent Integral. For further information please refer to the Fully-auto chemiluminescence immunoassay (CLIA) analyzer MAGLUMI Operating Instructions.

40µl	Sample, calibrator
------	--------------------

+80µl	ABEI Label
+20µl	Nano magnetic microbeads
15 min	Incubation
400µl	Cycle washing
3 s	Measurement

## DILUTION

Sample dilution by analyzer is not available in this reagent kit. Samples with concentrations above the measuring range can be diluted manually. After manual dilution, multiply the result by the dilution factor.

Please choose applicable diluents or ask SNIBE for advice before manual dilution must be processed.

## QUALITY CONTROL

- Observe quality control guidelines for medical laboratories
- Use suitable controls for in-house quality control. Controls should be run at least once every 24 hours when the test is in use, once per reagent kit and after every calibration. The control intervals should be adapted to each laboratory's individual requirements. Values obtained should fall within the defined ranges. Each laboratory should establish guidelines for corrective measures to be taken if values fall outside the range.

## LIMITATIONS OF THE PROCEDURE

### 1) Limitations

When basal LH values are interpreted, it must be taken into consideration that in healthy men, LH levels are subject to great spontaneous variations which may account for up to 20 pulses in 24 hours. Primary hypogonadism in men is not only associated with an elevation in mean LH concentrations but also with an increase in pulsatile secretion, while there are either no or only very sparse LH pulses, when the hypothalamus fails to secrete LHRH.

In women, a coincidentally measured mid-cycle gonadotropin peak may falsely point to primary ovarian insufficiency. A differentiation is possible when estradiol is determined at the same time. Therefore, a diagnosis should never be based on the determination of LH alone, but should always include other laboratory parameters and clinical findings.

### 2) Interfering Substances

No interference with test results is seen by concentrations of bilirubin < 66mg/dl, haemoglobin < 1000mg/dl or triglycerides < 1900mg/dl.

### 3) HAMA

Patient samples containing human anti-mouse antibodies (HAMA) may give falsely elevated or decreased values. Although HAMA-neutralising agents are added, extremely high HAMA serum concentrations may occasionally influence results.

### 4) High-Dose Hook

High dose hook is a phenomenon whereby very high level specimens may read within the dynamic range of the assay. For the MAGLUMI LH assay, No high-dose hook effect was seen for LH concentrations up to 3000 mIU/ml

## RESULTS

### 1) Calculation of Results

- The analyzer automatically calculates the LH concentration in each sample by means of a calibration curve which is generated by a 2-point calibration master curve procedure. The results are expressed in mIU/ml. For further information please refer to the Fully-auto chemiluminescence immunoassay (CLIA) analyzer MAGLUMI Operating Instructions.

### 2) Interpretation of Results

- Reference values:
 

Men:	1.1-25mIU/ml
Women:	follicular phase 1.2-12.5mIU/ml
	ovulatory period 12-82mIU/ml
	luteal phase 0.4-19mIU/ml

menopause 14-48mIU/ml

- Results may differ between laboratories due to variations in population and test method. If necessary, each laboratory should establish its own reference range.

## PERFORMANCE CHARACTERISTICS

### 1) Precision

Intra-assay coefficient of variation was evaluated on 3 different levels of control serum repeatedly measured 20 times in the same run, calculating the coefficient of variation.

Intra-assay precision			
Control	Mean(mIU/ml)	SD(mIU/ml)	CV%
Level 1	2.23	0.16	7.17%
Level 2	15.83	0.83	5.24%
Level 3	55.38	2.85	5.14%

Inter-assay coefficient of variation was evaluated on three batches of kits. Repeatedly measured 3 different levels of control serum 21 times, calculating the coefficient of variation.

Inter-assay precision			
Control	Mean(mIU/ml)	SD(mIU/ml)	CV%
Level 1	1.89	0.12	6.35%
Level 2	13.56	1.11	8.19%
Level 3	52.34	4.82	9.21%

### 2) Analytical Sensitivity

The sensitivity is defined as the concentration of LH equivalent to the mean RLU of 20 replicates of the zero standard plus two standard deviations corresponding to the concentration from the standard curve. The sensitivity is typically less than 0.5mIU/ml.

### 3) Specificity

The specificity of the LH assay system was assessed by measuring the apparent response of the assay to various potentially cross reactive analytes.

Compound	Concentration	Cross reactivity
FSH	150 mIU/ml	1.3%
HCG	500 mIU/ml	1%

### 4) Recovery

Consider calibrator high of known concentration as a sample, dilute it by 1:2 ratio with diluents, and measure its diluted concentration for 10 times. Then calculate the recovery of measured concentration and expected concentration. The recovery should be within 90% -110%.

Expected	Mean Measuring	Recovery
107.854 mIU/ml	112.684 mIU/ml	104%

### 5) Linearity

Use LH calibrator to prepare the six-point standard curve, measuring all points' RLU except point A, and then do four-parameter linear fitting in double logarithm coordinate, the absolute linear correlation coefficient(r) should be bigger than 0.9800.

Calibrator Point	Concentration mIU/ml	Absolute linear correlation coefficient (r)
A	0	
B	2	r=0.9940
C	5	
D	25	
E	75	
F	250	

### 6) Method comparison

A comparison of MAGLUMI LH(y) with a commercially available LH(x) using clinical samples gave the following correlations(mIU/ml):

Linear regression  
 $Y=0.96x+12.0$   
 $r=0.953$

Number of samples measured: 200

The sample concentrations were between 1.5 and 225.5 mIU/ml.

## REFERENCES

- Pagana, K. D. & Pagana, T. J. (© 2007). Mosby's Diagnostic and Laboratory Test Reference 8th Edition: Mosby, Inc., Saint Louis, MO. Pp 629-631.
- Wu, A. (© 2006). Tietz Clinical Guide to Laboratory Tests, 4th Edition: Saunders Elsevier, St. Louis, MO. Pp 694-697.
- Beastall GH, Ferguson KM, O'Reilly DSJ, Seth J, Sheridan B. Assays for follicle stimulating hormone and luteinizing hormone: Guidelines for the provision of a clinical biochemistry service. *Ann Clin Biochem.* 1987;24:246-262.
- Collip JB. William Henry Welch lectures: Some recent advances in physiology of anterior pituitary J. *MA Sinai Hosp.* 1934;1:28-71.
- Johnson MR, Carter G, Grint C, Lightmann SL. Relationship between ovarian steroids, gonadotropin and relating during the menstrual cycle. *Acta Endocrinol.* 1993;129/2:121-125.
- Runnebaum B, Rabe T. *Gynäkologische Endokrinologie und Fortpflanzungsmedizin* Springer Verlag 1994. Band 1:17,202-205,252-253, Band 2:350,360-362. ISBN 3-540-57345-3, ISBN 3-540-57347-x.
- Scott MG, Ladenson JH, Green ED, Gast MJ. Hormonal evaluation of female infertility and reproductive disorders. *Clin. Chem.* 1989;35:620-630.
- Van Ginkel LA, Loeber LG. Heterogeneity of human luteinizing hormone. Discrimination between acidic and basic preparations. *Acta Endocrinol.* 1989;121:73-82.
- Juhl UM, Rippegather G, Weller J, Zawta B. Important Facts on Reproduction Medicine/Fertility Diagnostics, Questions and Answers. 1994. Boehringer Mannheim, Cat.-No. 1322958.
- Tietz NW. *Clinical Guide To Laboratory Tests.* 3rd ed. Philadelphia, Pa: WB Saunders Co, 1995:410.
- Brzyski, R. and Jensen, J. (Revised 2007 March) *Female Reproductive Endocrinology, Introduction.* Merck Manual for Healthcare Professionals.
- Clarke, W. and Dufour, D. R., Editors (© 2006). *Contemporary Practice in Clinical Chemistry: AACC Press,* Washington, DC. Pp 360-361.
- Neely EK, Wilson DM, Lee PA, Stene M, Hintz RL (July 1995). Spontaneous serum gonadotropin concentrations in the evaluation of precocious puberty. *J Pediatric* 127(1):47-52 from PubMed.
- Genazzani AD, Petraglia F, Sgarbi L, Montanini V, Hartmann B, Surico N, Biolcati A, Volpe A, Genazzani AR. Difference of LH and FSH secretory characteristics and degree of concordance between postmenopausal and aging women. *Maturitas* 1997; 26: 133-138
- Huerta R, Malacara JM, Fajardo ME, Nava LE, Bocanegra A, Sanchez J. High-Frequency FSH and LH Pulses in Obese Menopausal Women. *Endocrine* 1997; 7 (3): 281-286
- Kim YK, Wasser SK, Fujimoto VY, Klein NA, Moore DE, Soules, MR. Utility of follicle stimulating hormone (FSH), luteinizing hormone (LH), oestradiol and FSH:LH ratio in predicting reproductive age in normal women. *Human Reproduction* 1997; 12 (6): 1152-1155
- Mahesh VB, Dhindsa DS, Anderson E, Kalra SP (Ed). *Regulation of Ovarian and Testicular Function.* Plenum Press, New York - London, 1988
- McCann SM, Kimura M, Walszewska A, Karanth S, Rettori V, Yu WH. Hypothalamic Control of FSH and LH by FSH-RF, LHRH, Cytokines, Leptin and Nitric Oxide. *Neuroimmunomodulation* 1998; 5: 193-202
- Pincus SM, Veldhuis JD, Mulligan H, Iranmanish A, Evans WS. Effects of age on the irregularity of LH and FSH serum concentrations in women and men. *Am J Physiol* 273 (Endocrinol Metab 36): E989-E995, 1997